

SINGLE USE SURGICAL PROFESSIONAL RANGE



CONTENTS

The Professional Range	4
Suctions	5-27
Suction Design	5
Suction Size Chart	6
Polished Suctions	7
Rounded Suctions	8
Fenestrated Suctions	9
Fenestrated & Rounded Suctions	10
Long Suctions	11-12
Teardrop Suctions	13
Fine Suctions	14
Soft Tip Suction	14
Retractor Suction	15
Dissector Suction	15
Micro Suction	16
Curved Micro Suction	16
Sinus Suctions	17
Nasal Aspirating Dissector	18
Schuknechts	18
Curved Suction	19
Poole Tip Suction	19
Poole Suction	19
Yankauer Suction	20
Sheffield Suction	20
Laser Suctions	21
Rounded Fenestrated Suction	22
Weir Suction	22
Insulated Suctions	23
Suction Tubing	23
Diathermy Suctions	24-26
Diathermy Leads	25
Sigmoidoscope Suction	27
Aspiration Tube	28
Aspiration Needle	28
Cannula	29-31
Tibbs Cannula	29
Cannula Straight	29
Cannula Curved	29
Angled Cannula	29

Sinus Trocar Cannula	30
Spackman Cannula	31
General Metal Instruments	32-38
Tischler-Morgan Biopsy Punch	32
Cushings Artery Forceps	33
Jefferson Retractor	33
Mollinsons Retractor	33
Hudson Brace	34
Gigli Saw	34
Yasargil Retractor	34
Fine Suture Forceps	35
Iris Scissors Straight & Curved	35
Kerrisons	36
Phlebectomy Instruments	37
Mosquito Forceps	37
Vein Hooks	37
Tracheal Dilator	38
Tilley Nasal Forceps	38
Thudichum Specula	38
Otology Instruments	39-42
Myringotomy Knife	39
Myringotomy Kit	39
Jobson Horne Probe	40
Cerumen Hook	40
Rosen Insertor	41
Micro Ear Forceps	41
Ear Specula	42
Pulse Lavage	43-44
Packaging	45-46
Index by Product Reference	47-48

SUCTIONS

CANNULAS

METAL INSTRUMENTS

OTOLOGY INSTRUMENTS

PULSE LAVAGE

September 2022

THE PROFESSIONAL RANGE



The problem: Debris found inside a reusable suction tube *after* cleaning and sterilization.

GUARANTEED CLEAN

All cannulated instruments are manufactured in the UK, using processes that ensure internal lumens are guaranteed clean and sterile to eliminate the risk of contamination from manufacturing debris.

Every instrument is delivered to you in individual sterile pouches for immediate use in the sterile field.

SUPERIOR QUALITY

Eakin Surgical strongly believes in value for money; it supplies the highest quality single use devices at competitive prices.

All products comply with industry standards as well as the exacting standards of professionals.

COST EFFECTIVE

When reusable instruments cannot be guaranteed clean, or are being lost, damaged and frequently replaced, single use instruments are a convenient cost effective option.



The solution: Single use suction tubes

CUSTOMER FOCUS

Meeting customer needs is at the heart of the company. Customer demand has driven the development of the extensive PROFESSIONAL range of suction tubes and other devices for many areas of surgery.

Eakin Surgical's friendly customer service team is always happy to work with you to find the right instrument to meet your requirements.

CONTACT US

Eakin Surgical Ltd
Greypoint
Cardiff Business Park
Cardiff
CF14 5WF
United Kingdom

Telephone +44 (0)2920 767 800
Facsimile +44 (0)2920 767 801
Email info@eakinsurgical.com
Website www.eakinsurgical.com

STANDARD SUCTION DESIGN

EXTENSIVE RANGE

The PROFESSIONAL range of suctions offers an extensive choice of length, lumen size and tip finish to enable you to find the right match for procedure, patient and surgeon.


Unless otherwise specified, all suctions have a 30° bend (a) to optimise visibility and this angle is fixed. The working length is measured from the bend to the tip (c).

ERGONOMIC HANDLE

A unique handle design provides optimum ergonomics with slim, smooth contours for comfort and improved visibility. Engineered to prevent blocking.

CONTROL HOLE



The suction control hole allows fine manipulation of suction with minimal noise. As standard the control hole is round but we also produce a range with a teardrop shaped control hole where more sensitive adjustment is needed (see page 13). There are also a small number of suctions that have no control hole in the handle, these will be indicated next to the listing.

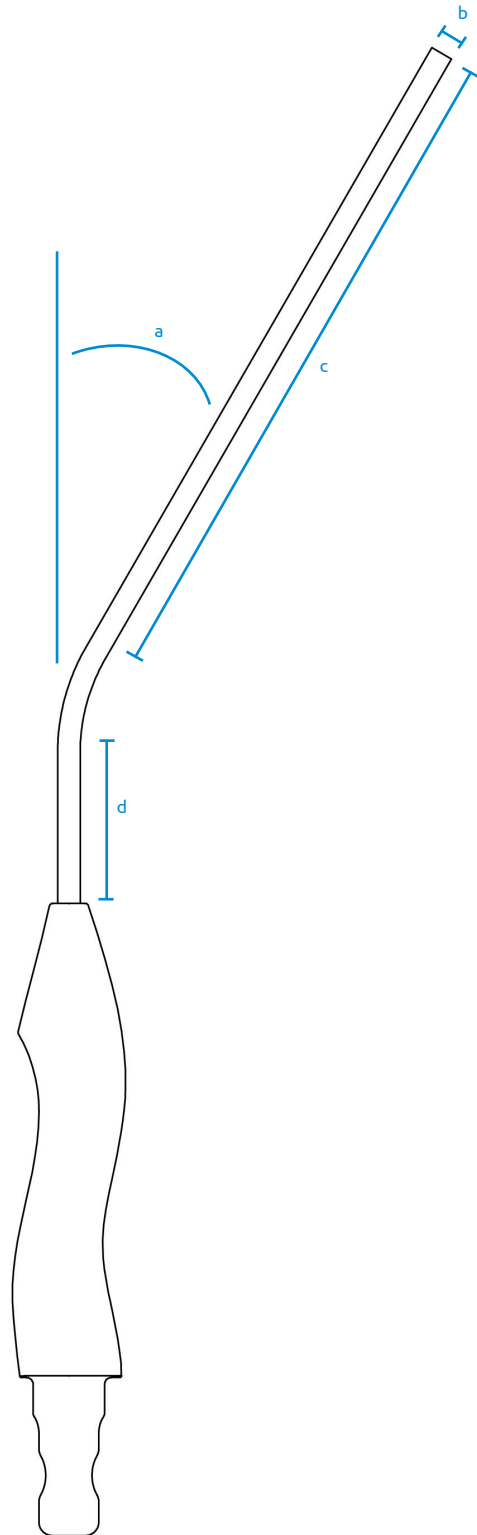

Ref.	Description
a	Bend angle
b	Lumen/Diameter
c	Working length
d	Length before the bend

Suction Size & Tip Style Chart

SUCTION

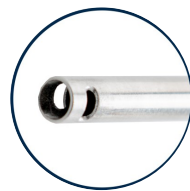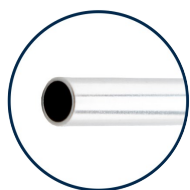

Tip	Polished								Rounded			Fenestrated			Rounded & Fenestrated		
fg	5	6	7	8	9	10	12	15	6	9	12	6	9	12	6	9	12
mm	1.7	2	2.3	2.7	3	3.3	4	5	2	3	4	2	3	4	2	3	4
5/6	5fgN100	6fgN100	7fgN100	8fgN100	9fgN100	10fgN100					S236	S086	S087	S160			
7									S102	S103							
9	5fgN140	6fgN140	7fgN140	8fgN140	9fgN140	10fgN140	S034	S340	S037	S014	S104	S002	S060		S091	S092	S093
13		S077			S076		S074	S390	S095	S096	S133	S055	S001	S042	S071	S094	S343
15					S005												
17							S199								S120	S075	
20		S031			S024		S021		S156	S319	S318				S168	S169	
23							S459										
27		S020			S011				S185	S182							
34							S460										
35					S036					S187	S284						
40							S109										
44							S461										
45		S258			S015					S183							
50							S132										
55					S111		S038				S184						
lumen	Working length (cm)																


Working length (cm)


Standard Range
Made to order
*Handle does not have a control hole

Polished Tips

All suctions have smooth, polished tips.

For rounded and fenestrated tips, please see pages 8-10. See also our range of specialist suctions on pages 12-27.

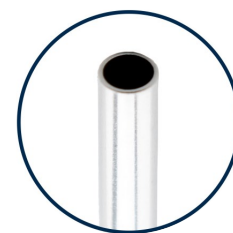

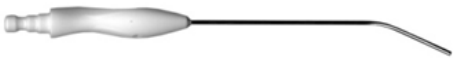	Description	Length	Lumen	Ref.
	Zoellner	5.5cm	5fg	5fgN100
			6fg	6fgN100
			7fg	7fgN100
	Mastoid	5.5cm	8fg	8fgN100
			9fg	9fgN100
			10fg	10fgN100
	Frazier	9cm	5fg	5fgN140
			6fg	6fgN140
			7fg	7fgN140
			8fg	8fgN140
			9fg	9fgN140
			10fg	10fgN140
			12fg	S034
			15fg	S340
		13cm	6fg	S077
			9fg	S076
			12fg	S074
			15fg	S390
		15cm	9fg	S005
		17cm	12fg	S199
	Magill	3.5cm	6fg	6fgF140
			9fg	9fgF140

Made to order

Rounded Tips

A range of suction tubes with a smooth rounded tip for atraumatic tissue contact.

For rounded tip suctions with a working length of 20cm and over, please see page 11.




	Description	Length	Lumen	Ref.
	Rounded Suction	5.5cm	12fg	S236
	Rounded Frazier	7cm	6fg	S102
			9fg	S103
	Rounded Frazier	9cm	6fg	S037
			9fg	S014
			12fg	S104
	Rounded Frazier	13cm	6fg	S095
			9fg	S096
			12fg	S133
	Rounded Magill	3.5cm	9fg	S046

Fenestrated Tips

A range of suctions with fenestrated tips to disperse aspiration and prevent suction trauma.

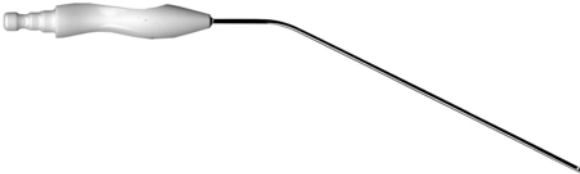

	Description	Length	Lumen	Ref.
	Fenestrated Suction	5cm	6fg	S086
			9fg	S087
			12fg	S160
	Fenestrated Suction	9cm	6fg	S002
			9fg	S060
	Fenestrated Suction	13cm	6fg	S055
			9fg	S001
			12fg	S042


Made to order

Fenestrated & Rounded Tips

Suctions that benefit from a tip that is both fenestrated and rounded.

For rounded, fenestrated suction tips with a working length of 20cm and over, please see page 12.

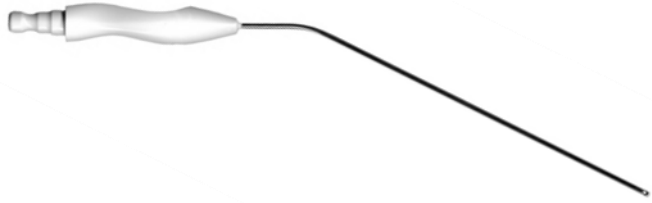


	Description	Length	Lumen	Ref.
	Fenestrated Rounded Frazier	5cm	6fg	S456
			9fg	S457
			12fg	S458
		9cm	6fg	S091
			9fg	S092
			12fg	S093
	Fenestrated Rounded Frazier	13cm	6fg	S071
			9fg	S094
			12fg	S343
	Fenestrated Rounded Suction	17cm	6fg	S120
			9fg	S075

Made to order

Long Suctions

A range of long suction tubes to optimize reach and visibility during ENT procedures.

- Brushed stainless steel prevents glare and reflection.
- Added benefit of rounded tips for extra atraumatic tissue contact.
- Strong high quality tube.

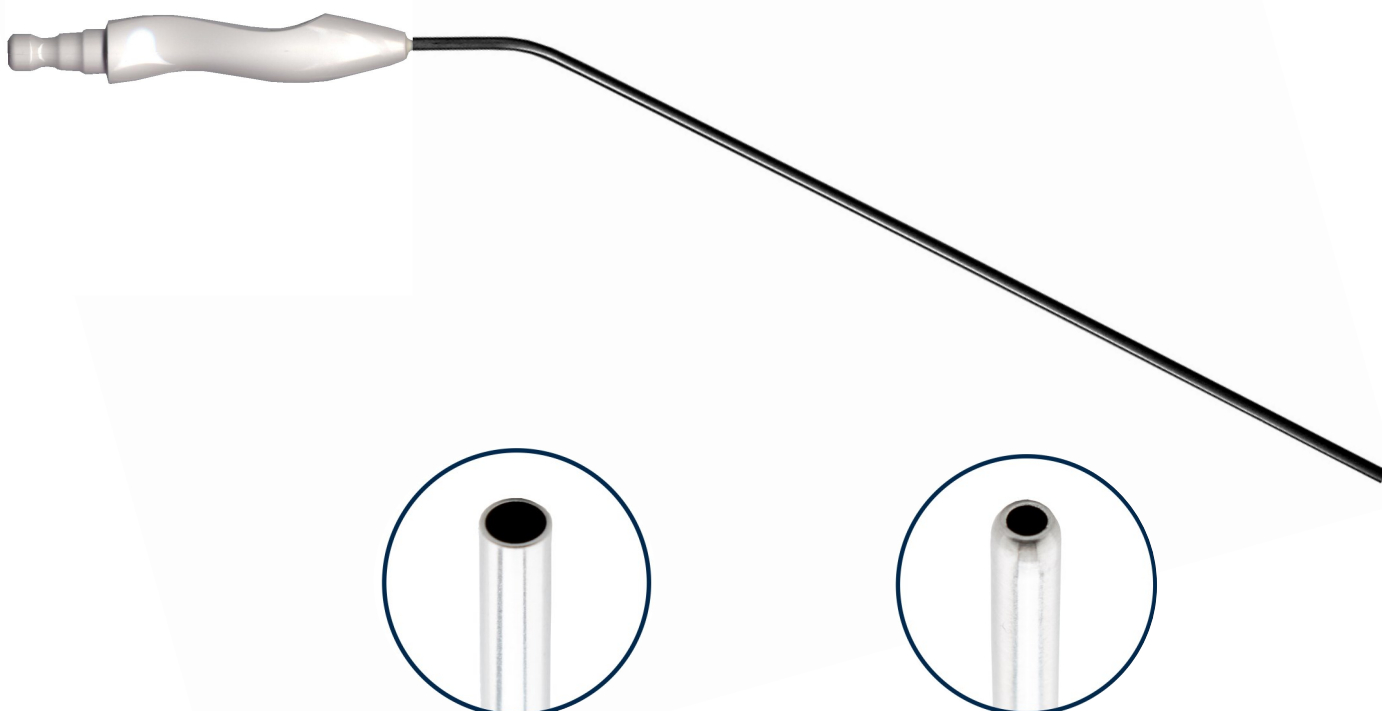


Working Length	Polished Tip			Rounded Tip		
	6fg ~ 2mm	9fg ~ 3mm	12fg ~ 4mm	6fg ~ 2mm	9fg ~ 3mm	12fg ~ 4mm
20 cm	S031	S024	S021	S156	S319	S318
23 cm			S459			
27 cm	S020	S011		S185	S182	
34 cm			S460			
35 cm		S036			S187	S284
40 cm			S109			
44 cm			S461			
45 cm	S258	S015			S183	
50 cm			S132			
55 cm		S111	S038			S184

Made to order

Long Suctions with 60° Bend

- 60° bend, provides a more acute angle than our standard suctions
- Brushed stainless steel prevents glare and reflection
- Strong high quality tube

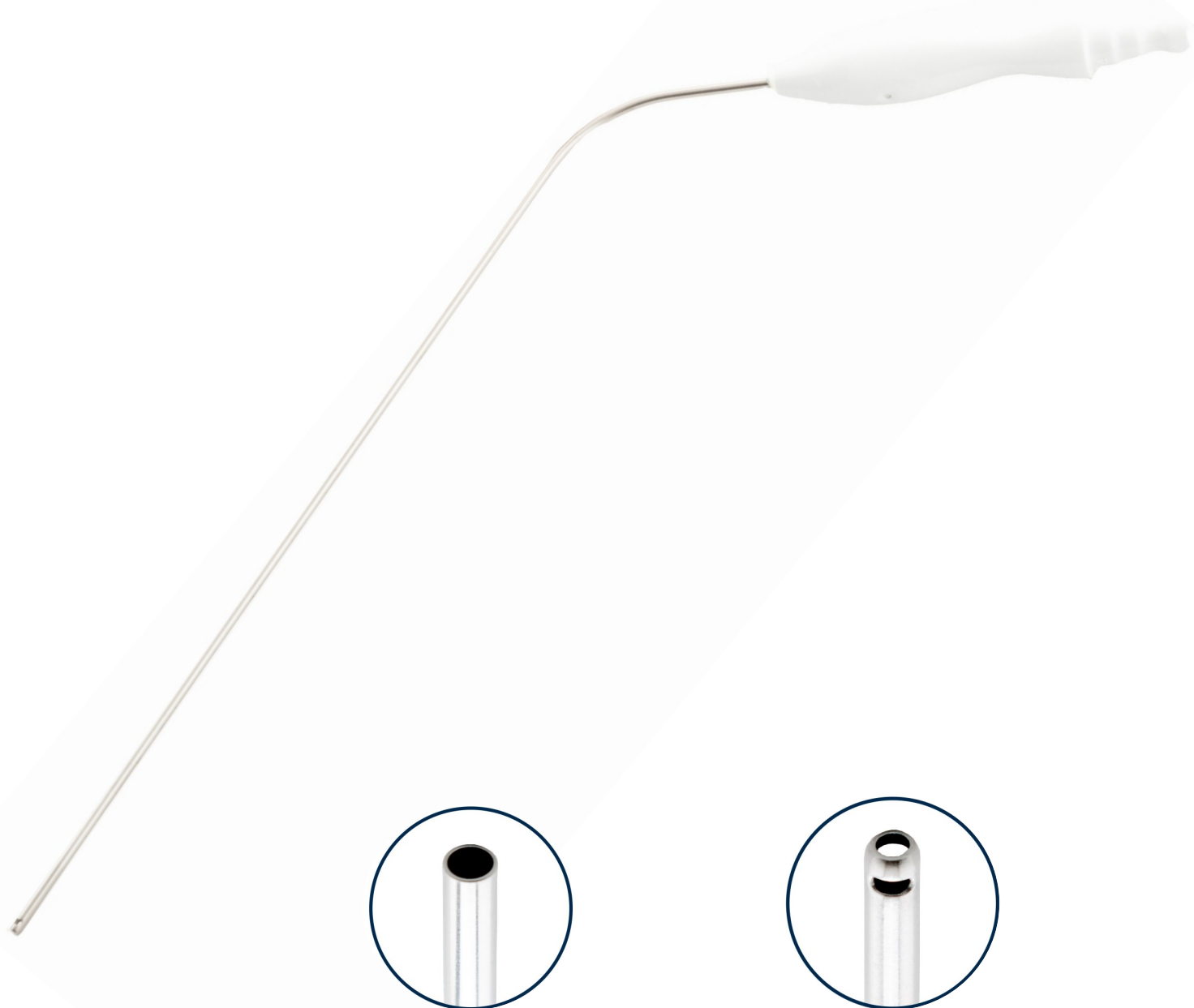


Working Length	Polished Tip	Rounded & Fenestrated Tip	
	12fg ~ 4mm	6fg ~ 2mm	9fg ~ 3mm
20 cm		S168	S169
30 cm	S131		

Made to order

Teardrop Control Suctions


- Teardrop shaped control hole for precise suction regulation
- Malleable tube gives the flexibility to adjust the angle by up to 50° for access to difficult to reach areas
- Available with a tapered tube to minimise blocking

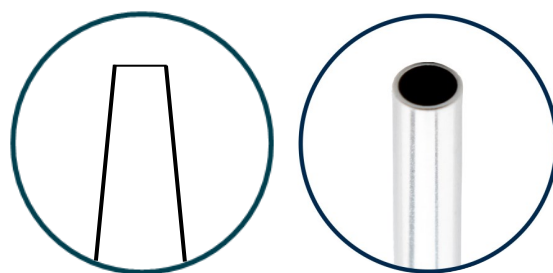

	Description	Length	Lumen	Ref.
	Teardrop Micro Suction	7cm	18g	S439
	Teardrop Suction Tapered	9cm	6fg	S190
	Teardrop Suction	9cm	7fg	S452
			8fg	S453
			9fg	S189
			10fg	S454
			12fg	S455
	Teardrop Suction	13cm	9fg	S192

Made to order

Fine Suctions

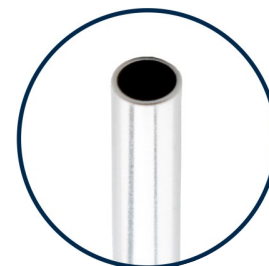

Reference	Description	Length	Lumen
S176	Suction 3fg 8cm	8cm	3fg
S177	Suction 5fg 8cm	8cm	5fg

[Made to order](#)

Soft Tip Suction

A suction with a soft silicone tip for very delicate procedures



Reference	Description	Length	Lumen
S321	Soft Tip 6fg 9cm Suction	9cm	6fg

Retractor Suction

A suction with a curved lip at its tip to allow simultaneous retraction and aspiration, available with teardrop control option.

The tube is flattened for easy access to the surgical site



Reference	Description	Length
S435	Retractor Suction 12fg 12cm 60°	12cm
S450	Teardrop Retractor Suction 12fg 12cm 60°	12cm

Made to order

Dissector Suction with Flattened Tip


A straight tube finished with a flattened polished tip




Reference	Description	Length
S436	Dissector Suction Flattened Tip 12fg	13cm

Micro Suction

The micro suction with an 18g fine tip

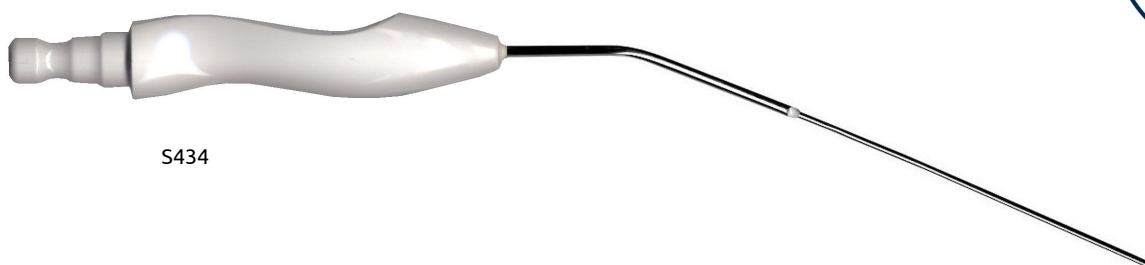


S434

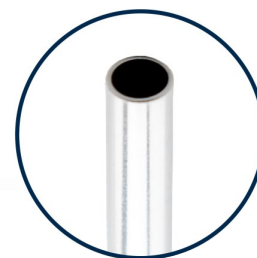

Curved Micro Suction

Anatomically designed to access harder to reach areas. Perfect for procedures performed through the nose, the smooth rounded tip allows for gentle, atraumatic dissection.

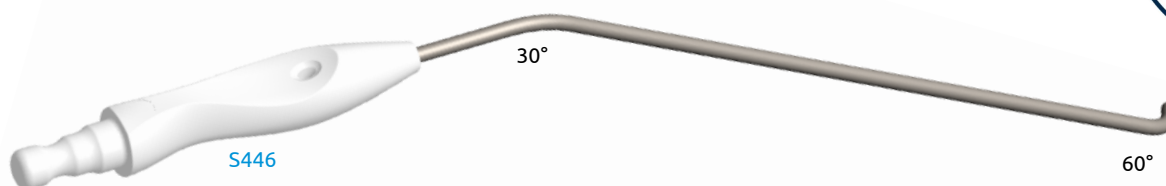

S446

60°



S447

Reference	Description	Length	Lumen
S434	Micro Suction 18g 7cm	7cm	18g
S446	Curved Micro Suction 9fg 15cm 30° Bend Rounded Tip 60° Curve (Left)	15cm	9fg
S447	Curved Micro Suction 9fg 15cm 30° Bend Rounded Tip 60° Curve (Right)	15cm	9fg

Made to order

Sinus Suctions

A range of suctions with curved tips for maximum visibility and access in difficult to reach areas



Extended Range



Gentle bend for improved access

With control hole only

60° Curve Ball End

9fg—3mm

S444



	Curve Rounded Tip		Curve Ball End		Curve Rounded Tip		Curve Ball End	
	70°		70°	65°	90°		90°	
	9fg ~ 3mm	12fg ~ 4mm	9fg ~ 3mm	12fg ~ 4mm	9fg ~ 3mm	12fg ~ 4mm	9fg ~ 3mm	12fg ~ 4mm
	<i>Upward Curve</i>		S467	S468				
<i>With Control Hole</i>	S012	S013	S228	S229	S040	S041	S230	S231
<i>No Control Hole</i>		S163						

Made to order

Nasal Aspirating Dissector

The flattened curved tip of this instrument makes it ideal for access and elevation in the nasal cavity. A small fenestration on both sides of the tip allows simultaneous suction



Reference	Description	Length
S146	Nasal Aspirating Dissector	12cm

Schuknechts

Schuknechts provide fine end aspiration



Reference	Description	Length	Lumen
S007	Schuknecht 20g	6cm	20g
S008	Schuknecht 24g	6cm	24g
S072	Schuknecht 18g	6cm	18g
S157	Long Schuknecht 20g	13cm	20g

Made to order

Curved Suction

This suction has a smooth curvature with rounded and fenestrated tip for atraumatic and dispersed aspiration



Reference	Description	Length	Lumen
S026	Curved Suction 12fg	14cm	12fg

Poole Tip Suction



Reference	Description	Length	Lumen
S440	Poole Tip Suction 12fg 13cm Rounded 4x4x2mm Fenestrations	13cm	12fg

Poole Suction

The Poole suction has multiple holes around the sheath for dispersed suction and reduced risk of blockage. It also offers the option for finer suction with a simple bayonet twist action




Reference	Description	Length	Lumen
S010	Poole Suction	26cm (when closed)	5mm

Yankauer Suction

New

Available in ball and fenestrated tip style variants and ensures reliability throughout the entire procedure

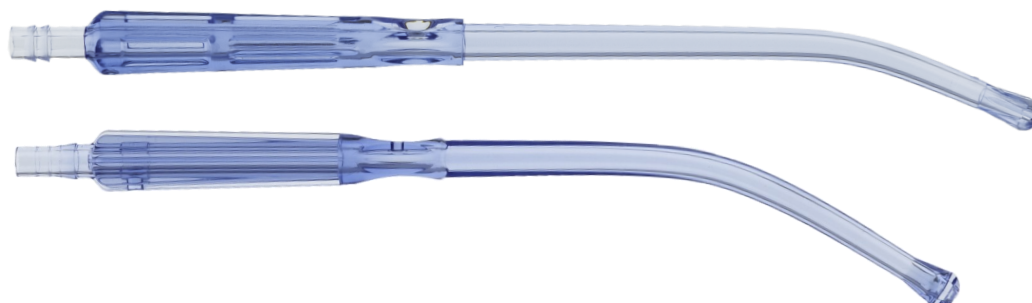


Reference	Description
GC1233K3W	Yankauer Suction Ball Tip
GC1215A3W	Yankauer Suction Fenestrated Tip

Yankauer Suctions are under distribution agreement of GCMEDICA Enterprise Ltd.

Sheffield Suction

Suction with novel clip to allow attachment to a laryngoscope during endolaryngeal laser surgery.

Enables effective smoke plume evacuation and is detachable in the instance of a bleed

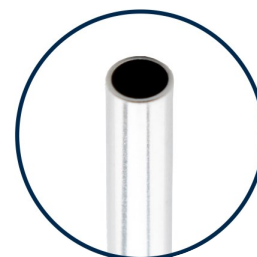

“ We routinely use these sucker clips to perform all endolaryngeal [laser] surgery and overall, feel they are simple, easy to use and extremely efficient in suctioning fumes. ”

Mr Jabreel, Consultant ENT Surgeon


Mr Manickavasagam, Consultant ENT Surgeon

Reference	Description	Length	Lumen
S240	Sheffield Suction 9fg 15cm 90° Bend with Clip	15cm	9fg

Laser Suctions

Laser suctions combine two instruments in one; a 0.9mm fibre optic laser port bonded to the suction tube allows simultaneous aspiration with cutting and coagulation

For use with 600 micron core laser fibre with 850 micron sheath

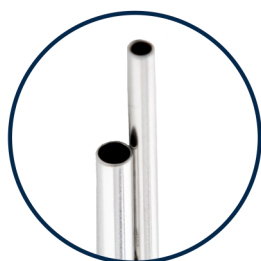

Extended Laser Channel



Flush Finish Laser Channel



Extended Angulated Laser Channel

	Description	Length	Lumen	Ref.
	Laser Suction 6fg 9cm with Extended Laser Channel	9cm	6fg	S110
	Laser Suction with Flush Finish	9cm	6fg	S165
		20cm	6fg	S139
	Fairley Laser Suction 12cm with Extended Angulated Tip	12cm	6fg	S285

Made to order

Rounded Suction with 4 Fenestrations

This ultra delicate suction has 4 fenestrations that are 1mm in diameter at the tip



Reference	Description	Length	Lumen
S147	Rounded Suction 12fg 9cm 4x1mm Fenestrations	9cm	12fg
S445	Rounded Suction 12fg 30cm 4x1mm Fenestrations 60°	30cm	12fg

Weir Suction

A suction with a flattened, semi-sharp open end to combine minor dissection with aspiration




Reference	Description	Length
S009	Weir Suction 12fg	13cm

Made to order

Insulated Suctions

A range of suctions with protective insulation to allow use alongside diathermy

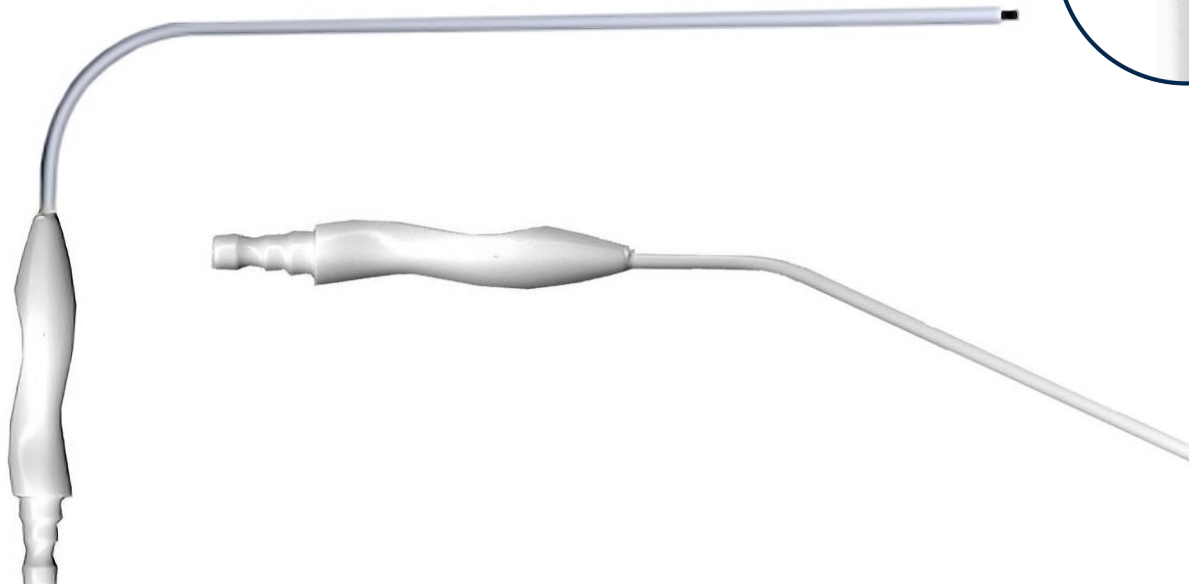

Reference	Description	Length	Lumen
S432	Insulated Suction 9fg 17cm 90°	17cm	9fg
S433	Insulated Suction 6fg 17cm 90°	17cm	6fg
S025	Insulated Frazier 9fg 9cm	9cm	9fg

Made to order

Flexible Connection Tubing

Made from silicone for flexibility and resistance to kinking



Connector

Reference	Description	Includes Connector
S066	Flexible Connection Tubing 8mm 50cm	No
S099	Flexible Connection Tubing 8mm 300cm	Yes
S101	Flexible Connection Tubing 8mm 150cm	Yes

Diathermy Suctions

As a 2 in 1 instrument, diathermy suctions create a clearer view and a less invasive procedure by minimising the number of instruments in the field of vision.

The suction helps remove harmful smoke from theatre. The 4mm pin connects to diathermy lead (see opposite page for reusable leads).



Features on-tip smoke evacuation as recommended by HSE



Working Length	Polished Tip		Rounded Tip	Fenestrated Tip	Rounded & Fenestrated Tip
	9fg ~3mm	12fg ~4mm	9fg ~ 3mm	9fg ~3mm	9fg ~ 3mm
9 cm	S117		S136		
13 cm	S116				S119
19 cm		S195			
20 cm	S115				
22 cm				S113	
27 cm	S114				

Diathermy Suctions

Featuring a 70° bend and a standard polished tip.



Diathermy Leads

Competitively priced reusable leads for use with Single Use Surgical 'Professional' Diathermy suctions



Reference	Description	Length	Lumen
S438	Diathermy Suction 9fg 17cm 70°	17cm	9fg
1011-35	Diathermy lead 4mm pin 4mm socket	3m	n/a
1012-35	Diathermy lead 8mm pin 4mm socket	3m	n/a

Made to order

Diathermy Suctions

Anatomically designed to access harder to reach areas often associated with Submucosal Diathermy (SMD).

Combined diathermy and suction for targeted and effective smoke extraction

- Ø3mm—6fg 10cm
- 40° bend
- 4mm connection pin
- Latex & DEHP free



Built in smoke extraction



Features on-tip smoke evacuation as recommended by HSE


“ I have used the Diathermy Abbey Needle with Suction, and I like it. It is quicker and less traumatic to use. It feels more robust than the usual needle, so you can push the turbinate laterally as well, creating a bit more airway and avoiding accidental diathermy of the septum. This reduces the chance of complications such as intra-nasal adhesion occurring. Simultaneous clearance of the smoke means you don't need to alternate between the diathermy needle and the Zoellner sucker, which can often cause minor bleeding. It allows good vision throughout the procedure. ”

David Hill, Consultant ENT Surgeon

Reference	Description	Length	Lumen
S448	Diathermy Abbey Needle with Suction	10cm	6fg

Sigmoidoscope Suction

A rigid sigmoidoscope suction without a suction control hole

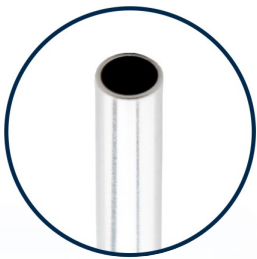

No fine control hole in handle to prevent leakage



Reference	Description	Length	Lumen
S043	Sigmoidoscope Suction 12fg 44cm	44cm	12fg

Made to order

Aspiration Needles

Aspiration needles with a sharp bevelled tip for safe and controlled penetration when aspirating fluids from cysts and organs. With Luer Lock.


Please note, product colour may vary from that pictured.

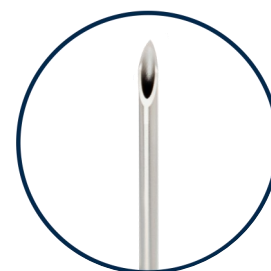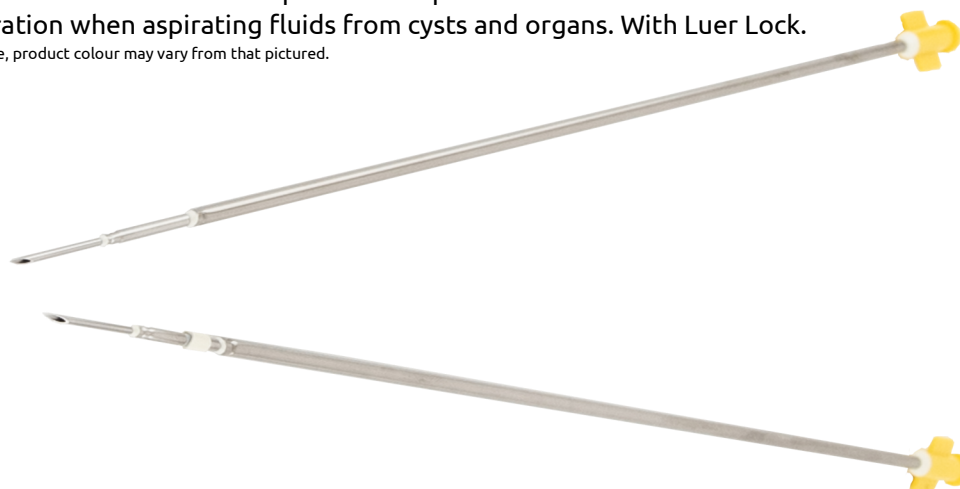

Reference	Description	Length
S019	Aspiration Needle 3mm 32cm	32cm
S080	Aspiration Needle 5mm 32cm	32cm

Aspiration Tube


Atraumatic blunt ended tube with Luer Lock

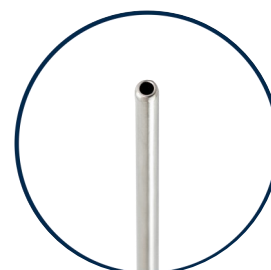

Reference	Description	Length	Lumen
S112	Blunt Aspiration Tube 4mm 33cm	33cm	4mm

Tibbs Cannula

The Tibbs hybrid cannula combines traditional olive and cone tips

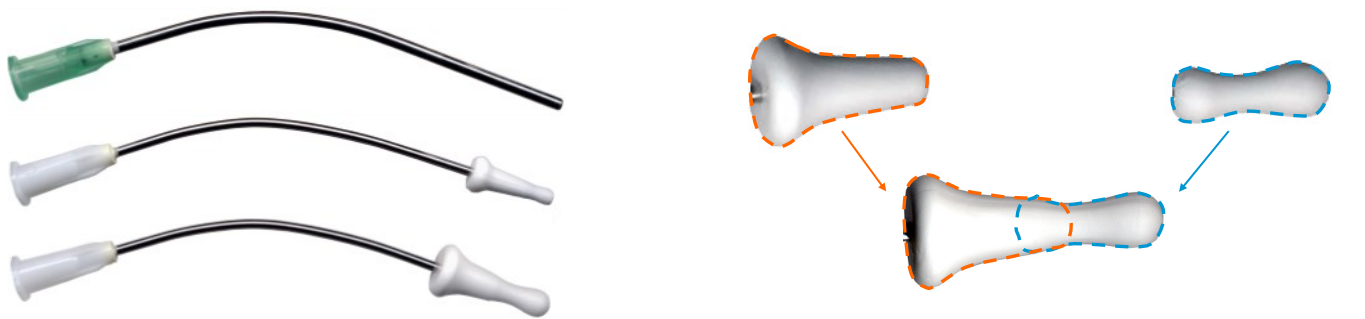

Reference	Description	Length
S044	Tibbs Cannula Set of 3	7cm

Straight, Curved & Angled Cannulas




Reference	Description	Length	Lumen
S214	6fg Cannula Straight	7cm	6fg
S213	6fg Cannula Curved Set of 2	7cm	6fg
S028	Angled Cannula 19g	2.5cm	19g

Sinus Trocar Cannula

The Sinus Trocar and Cannula - also known as an antrum cannula - consists of a removable cannula over a 11.2cm trocar to allow washout of the sinuses.

Trocar 16g 11.2cm
Cannula 14g 10.7cm with Luer Lock

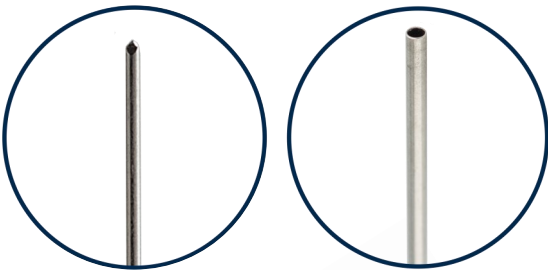

Reference	Description
S081	Sinus Trocar Cannula

Spackman Cannula

The single use Spackman Cannula combats concerns over cleanliness with the added benefit of an easy to use self-locking clamp. Available with a fixed or repositionable acorn.


For hydrotubation dye tests and uterine manipulation during laparoscopy.



Fixed Acorn



Moveable Acorn

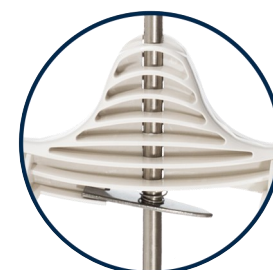Rounded &
Fenestrated Tip

Self-Locking Clamp



Reference	Description	Length	Lumen
S029	Spackman Cannula	32cm	12fg
S030	Spackman Cannula Moving Acorn	32cm	12fg

Tischler-Morgan Biopsy Punches

The biopsy punches come with either a straight tip (pictured right) or a tip that is angled down (pictured below) to allow surgical staff to select the best option for their specific patient and procedure.

- Single use eliminates patient cross infection risk, biopsy sample cross contamination and guarantees quality
- Tip locking system safely secures sample
- Ergonomic hand-polished edges
- Quality & functionality matches reusables
- Bite: 7mm x 3.5mm
- Working Length: 200mm



Reference	Description
S263	Tischler-Morgan Biopsy Punch 200mm, Straight Tip
S320	Tischler-Morgan Biopsy Punch 200mm, Tip Down


Cushings Artery Forceps



Reference	Description
S205	Cushings Artery Forceps Curved, Half Serrated Jaws, 150mm

Made to order

Jefferson Retractor

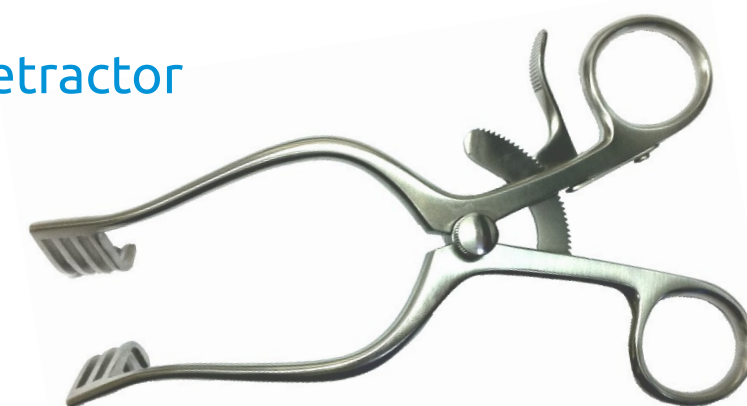


Reference	Description
S206	Jefferson Retractor, Self-Retaining, Curved, Blunt Pointed Blades, 125mm

Mollinsons Retractor



Reference	Description
S207	Mollinsons Retractor, Self-Retaining, Curved, 4/4 Teeth, Adult, 165mm

Hudson Brace, Gigli Saw and Yasargil Retractor

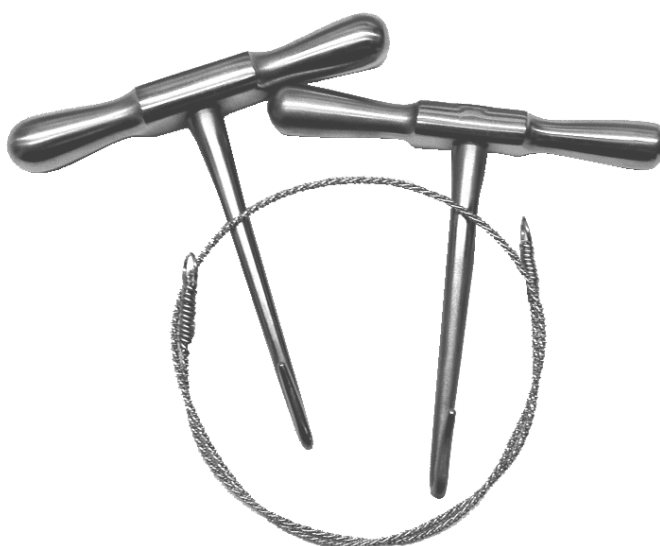

Reference	Description
S211	Hudson Brace, Standard Size with Universal Fitting, 265mm
S212	Gigli Saw 300mm and Handles, Pair, 55mm with Guide
S419	Yasargil Retractor (not pictured)

Made to order

Fine Suture Forceps



Reference	Description
S257	Fine Suture Forceps

Made to order

Iris Scissors Straight & Curved



Reference	Description
S255	Iris Scissors Straight
S256	Iris Scissors Curved

Kerrisons

- Designed for single patient use to eliminate cross infection and guarantee a sharp cut
- Ergonomic hand-polished edges
- Quality & functionality matches reusables



Bite	45° Degree Upward Cutting (Delicate Use)	45° Degree Upward Cutting	90° Degree
1mm	S161	S125	
2mm	S261	S126	S367
3mm	S262	S127	
4mm		S137	
5mm		S138	

Made to order

Vein Hooks

A cost effective alternative to reusable equivalents. Made from high grade stainless steel they have a strong Oesch style hook which can withstand rigorous handling, yet is fine enough for insertion through small incisions



“ We have recently started using phlebectomy hooks and mosquito clips from Eakin Surgical Ltd at all of our clinics. The first thing we noticed was the very good quality. For disposable surgical instruments, the feel and performance of the hooks and the clips is excellent. As the quality of our work depends upon the quality of the instruments that we use, we are delighted to have found Eakin Surgical Ltd. ”

Professor Mark Whiteley, Venous Surgeon and owner of the Whiteley Clinic

Reference	Description
S090	Medium Right Hand Vein Hook (No. 2)
S106	Small Right Hand Vein Hook (No. 3)
S134	Medium Left Hand Vein Hook (No. 2)
S135	Small Left Hand Vein Hook (No. 3)

Mosquito Forceps

Made to order

Single use mosquito forceps provide convenience and consistent quality



Reference	Description
S098	Mosquito Forceps Curved Set of 2

Made to order

Tracheal Dilator



Reference	Description
S181	Tracheal Dilator

Made to order

Tilley Nasal Forceps



Reference	Description
S180	Tilley Nasal Forceps Short
S380	Tilley Nasal Forceps Long

Thudichum Specula



Reference	Description
S179	Thudichum Specula Size 5
S238	Thudichum Specula Size 4
S239	Thudichum Specula Size 3

Made to order

Myringotomy Kit

These all in one procedure packs contain sterilised, easy to use disposable instruments that are packed in a sterile peel pouch to save on cleaning and assembly costs



Reference	Description
S301	Myringotomy Kit 1 (Micro Forcecp, Knife, Jobson Horne Probe, Rosen Inserter, Ear Specula 3/4/5/6mm, Zoellner Suction Tube, 18g Zoellner Tip, 22g Zoellner Tip)

Myringotomy Knife

Lancet used to pierce the tympanic membrane. Bendable shaft for ease, versatility and to prevent visual obstruction. Handle is shaped to aid orientation of the blade.




Reference	Description	Total length	Handle length
S310	Myringotomy Knife	157mm	95mm

Jobson Horne Probe

Designed to allow gentle and easy extraction of wax or foreign bodies from the ear canal. The probe has a Curette loop at one end and a threaded section at the other that can be used to hold wrapped cotton wool.


The handle section in the centre features flattened sides to aid grip and prevent unintentional twisting. Made from carbon fibre filled nylon for strength and flexibility.

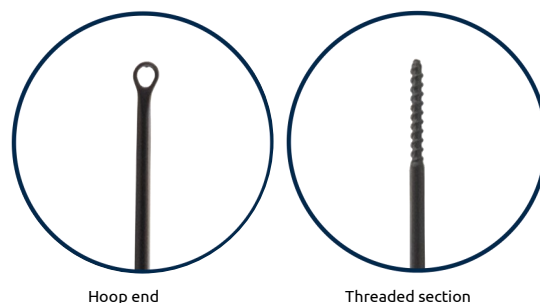

Reference	Description	Length	Loop
S303	Jobson Horne Probe	150mm	3mm

Cerumen Hook

Disposable cerumen hook for extracting hard wax and other obstacles from the ear canal. The hook shape has been designed to offer good leverage whilst the flexible carbon fibre helps reduce trauma to the ear.

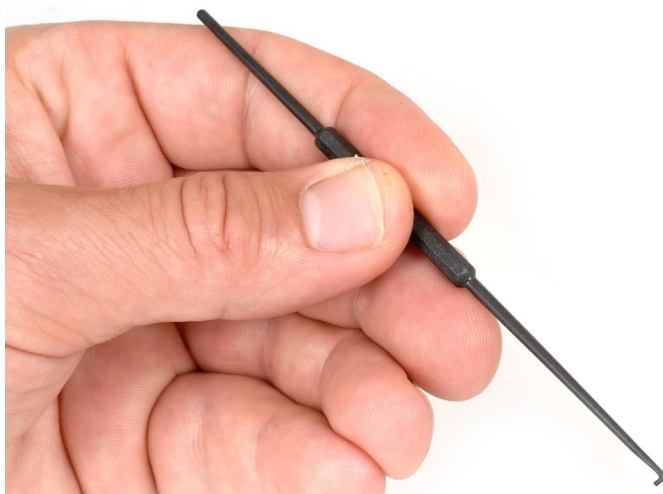


Reference	Description	Length	Hook Length
S307	Cerumen Hook	135mm	2mm

Made to order

Rosen Inserter

Developed with surgeons to aid grommet insertion. The geometry also proves favourable for removal of objects in the ear.

Manufactured from carbon fibre filled nylon for strength and flexibility.

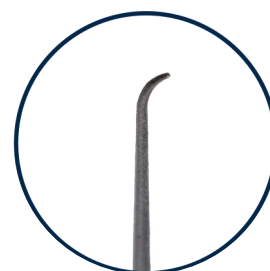

Reference	Description	Length
S308	Rosen Inserter	143mm

Made to order

Micro Ear Forceps

High quality micro ear forceps—also known as crocodile forceps—that effectively remove foreign bodies or are used to insert grommets. Feature serrated jaws to easily grip particles and fine shaft for access to the inner ear.



Reference	Description	Shaft length
S221	Micro Ear Forceps 5mm jaws	85mm

Ear Specula

Advanced design provides maximum patient comfort and superior operational vision

Matt finish and internal ridges reduce glare and reflection from the microscope light onto the internal surface of the specula. Whilst the end is slightly oval to conform to the ear canal

Soft and flexible material to aid patient comfort

Serrated lateral end aids finger tip grip



Reference	Description
S304	Ear Specula 4mm
S305	Ear Specula 5mm
S306	Ear Specula 6mm
S309	Ear Specula 3mm

[Made to order](#)

Pulse Lavage Kit

A disposable pulse lavage system with a built-in battery, for controlled pulsatile irrigation and suction

- Built-in battery pack - easily removed for safe disposal**, cutting the weight and cost of clinical waste
- Long coaxial tip provides extra reach
- 30mm cone tip
- Ergonomic handle

Technical Specification

- High speed flow rate up to 1,000ml/M
- Normal speed flow rate up to 700ml/M
- Pressure up to 15 PSI



Includes 2 additional tips



Built-in battery pack

Quick-locking clip
securely holds tips


Reference	Description
S353	Pulse Lavage Kit 1 includes Pulse Lavage, 30mm Cone Tip, Long Coaxial Canal (Femoral) Tip
Accessories - see next page for images	
BT 04 Z	Coaxial Canal (Femoral) Brush Tip
CT 02 Z	Coaxial Canal (Femoral) Tip
CTL 02 Z	Long Coaxial Canal (Femoral) Tip
LT 01 Z	75mm Cone Tip
MT 01 Z	50mm Cone Tip
ST 01 Z	30mm Cone Tip
RT23012Z	Hip Revision Tip

Made to order

** Ensure batteries are disposed of safely and in accordance with all current national and international waste disposal directives

Pulse Lavage - Accessories

CT 02 Z Coaxial
Canal (Femoral)
Tip

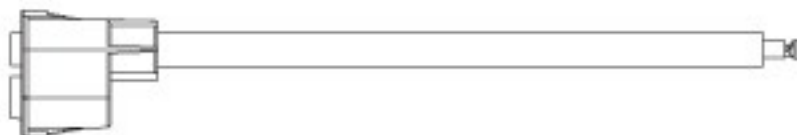

CTL 02 Z Long
Coaxial Canal
(Femoral) Tip




BT 04 Z Coaxial
Canal (Femoral)
Brush Tip




ST 01 Z 30mm
Cone Tip




MT 01 Z 50mm
Cone Tip

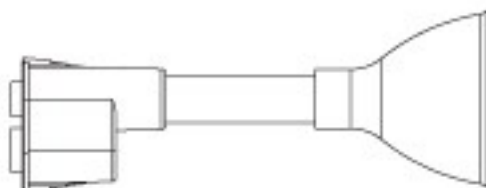

LT 01 Z 75mm
Cone Tip

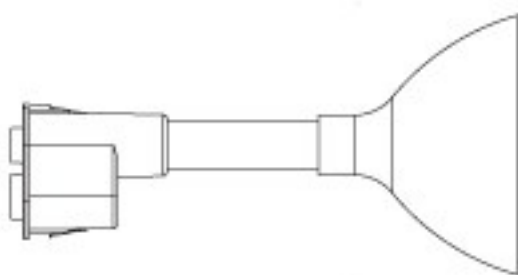

RT23012Z Hip
Revision tip

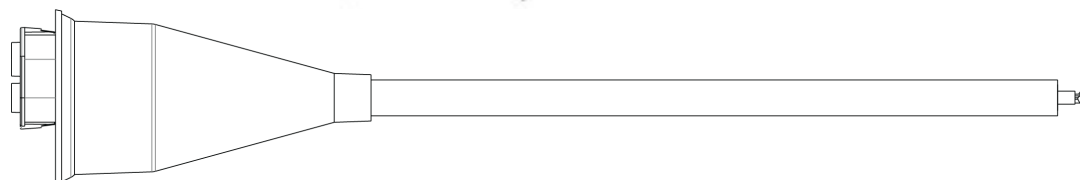

Pulse Lavage Tip accessories are under distribution agreement of Dongguan Kaiser Technology Co Ltd.

Made to order

Packaging

STERILE POUCHES

All sterile instruments are supplied in individual peelable pouches to protect them during storage until needed.

DISPENSING BOXES

Each PROFESSIONAL box has a removable lid for convenient dispensing of sterile pouches even when boxes are stacked. The front label on the boxes has a barcode for quick identification and stock management.

CARE FOR THE ENVIRONMENT

Eakin Surgical strive to protect the environment through responsible and sustainable operations. All cardboard packaging is recyclable. For information on recycling of metal instruments please ask your representative

Packaging Symbol Key

 Manufacturer	 Sterilized using ethylene oxide	 Non-sterile
 Catalogue number	 Do not re-sterilize	 Keep dry
 Batch code	 Single barrier system with protective packaging outside	 Unique Device Identifier
 Do not use if package is damaged and consult instructions for use	 Consult instructions for use or consult electronic instructions for use	 Latex Free
 Do not re-use	 Medical Device	 Temperature limit
 Caution	 Use-by-date	 Keep away from sunlight
 Date of Manufacture	 Authorized representative in the European Community/European Union	 Refer to instruction manual/booklet
 Quantity of content	 Single sterile barrier system	 Protect from heat and radioactive sources
 Indicates a recovery/recyclable package or packaging material	 Alkaline Battery (recyclable)	Recycle in accordance with The Waste Electrical and Electronic Equipment Directive (WEEE Directive) European Community directive 2012/19/EU
 CAUTION: US Federal law restricts this device to sale by or on the order of the physician		

 CE mark including SGS (notified body) identification number. Product conforms with European Requirements

 UKCA mark including SGS (notified body) identification number. Product conforms with UK Requirements

INDEX BY PRODUCT REFERENCE

Reference	Description	Page
1011-35	3m Diathermy lead 4mm pin, 4mm socket.	25
1012-35	3m Diathermy lead, 8mm pin, 4mm socket	25
10fgN100	Mastoid Suction 10fg	7
10fgN140	Frazier Suction 10fg 9cm	7
5fgN100	Zoellner Suction 5fg	7
5fgN140	Fine Frazier Suction 5fg 9cm	7
6fgF140	Fine Magill 6fg	7
6fgN100	Zoellner 6fg	7
6fgN140	Fine Frazier 6fg 9cm	7
7fgN100	Zoellner Suction 7fg	7
7fgN140	Fine Frazier Suction 7fg 9cm	7
8fgN100	Mastoid Suction 8fg	7
8fgN140	Frazier Suction 8fg 9cm	7
9fgF140	Magill Suction 9fg	7
9fgN100	Mastoid Suction 9fg	7
9fgN140	Frazier Suction 9fg 9cm	7
BT 04 Z	Coaxial Canal (Femoral) Brush Tip	44
CT 02 Z	Coaxial Canal (Femoral) Tip	44
CTL 02 Z	Long Coaxial Canal (Femoral) Tip	44
GC1215A3W	Yankauer Suction Ball Tip	20
GC1233K3W	Yankauer Suction Fenestrated Tip	20
LT 01 Z	75mm Cone Tip	44
MT 01 Z	50mm Cone Tip	44
RT23012Z	Hip Revision Tip	44
ST 01 Z	30mm Cone Tip	44
S001	Fenestrated Suction 9fg 13cm	9
S002	Fenestrated Suction 6fg 9cm	9
S005	Suction 9fg 15cm	7
S007	Schuknecht 20g 6cm	18
S008	Schuknecht 24g 6cm	18
S009	Weir Suction 12fg	22
S010	Poole Suction	19
S011	Long Suction 9fg 27cm	11
S012	Sinus 9fg	17
S013	Sinus 12fg	17
S014	Rounded Frazier 9fg 9cm	8
S015	Long Suction 9fg 45cm	11
S019	Aspiration Needle 3mm 32cm	28
S020	Long Suction 6fg 27cm	11
S021	Long Suction 12fg 20cm	11
S024	Long Suction 9fg 20cm	11
S025	Insulated Frazier 9fg 9cm	23
S026	Curved Suction 12fg	19
S028	Angled Cannula 19g	29
S029	Spackman Cannula	31
S030	Spackman Cannula Moving Acorn	31
S031	Long Suction 6fg 20cm	11
S034	Large Frazier 12fg 9cm	7
S036	Long Suction 9fg 35cm	11
S037	Rounded Fine Frazier 6fg 9cm	8
S038	Long Suction 12fg 55cm	11
S040	Sinus 9fg Small Bend	17
S041	Sinus 12fg Small Bend	17
S042	Fenestrated Suction 12fg 13cm	9
S043	Sigmoidoscope Suction 12fg 44cm	27
S044	Tibbs Cannula Set of 3	29
S046	Rounded Magill 9fg	8

Reference	Description	Page
S055	Fenestrated Suction 6fg 13cm	9
S060	Fenestrated Frazier 9fg 9cm	9
S066	Flexible Connection Tubing 8mm 50cm	23
S071	Fenestrated Rounded Frazier 6fg 13cm	10
S072	Schuknecht 18g 6cm	18
S074	Suction 12fg 13cm	7
S075	Fenestrated Rounded Suction 9fg 17cm	10
S076	Suction 9fg 13cm	7
S077	Suction 6fg 13cm	7
S080	Aspiration Needle 5mm 32cm	28
S081	Sinus Trocar/Cannula	30
S086	Fenestrated Suction 6fg 5cm	9
S087	Fenestrated Suction 9fg 5cm	9
S090	Medium Right Hand Vein Hook (No. 2)	37
S091	Fenestrated Rounded Frazier 6fg 9cm	10
S092	Fenestrated Rounded Frazier 9fg 9cm	10
S093	Fenestrated Rounded Frazier 12fg 9cm	10
S094	Fenestrated Rounded Frazier 9fg 13cm	10
S095	Rounded Frazier 6fg 13cm	8
S096	Rounded Frazier 9fg 13cm	8
S098	Mosquito Forceps Curved Set of 2	37
S099	Flexible Connection Tubing 8mm 300cm	23
S101	Flexible Connection Tubing 8mm 150cm	23
S102	Rounded Frazier 6fg 7cm	8
S103	Rounded Frazier 9fg 7cm	8
S104	Rounded Frazier 12fg 9cm	8
S106	Small Right Hand Vein Hook (No. 3)	37
S109	Long Suction 12fg 40cm	11
S110	Laser Suction 6fg 9cm	21
S111	Long Suction 9fg 55cm	11
S112	Blunt Aspiration Tube 4mm 33cm	28
S113	Diathermy Suction 9fg 22cm Fenestrated	24
S114	Diathermy Suction 9fg 27cm	24
S115	Diathermy Suction 9fg 20cm	24
S116	Diathermy Suction 9fg 13cm	24
S117	Diathermy Suction 9fg 9cm Frazier	24
S119	Diathermy Suction 9fg 13cm Fenestrated Rounded	24
S120	Fenestrated Rounded Suction 6fg 17cm	10
S125	Kerrison Punch 1mm x 200mm @45°	36
S126	Kerrison Punch 2mm x 200mm @45°	36
S127	Kerrison Punch 3mm x 200mm @45°	36
S131	Long Suction 12fg 30cm @60°	12
S132	Long Suction 12fg 50cm	11
S133	Rounded Suction 12fg 13cm	8
S134	Medium Left Hand Vein Hook (No. 2)	37
S135	Small Left Hand Vein Hook (No. 3)	37
S136	Diathermy Suction 9fg 9cm Rounded	24
S137	Kerrison Punch 4mm x 200mm @45°	36
S138	Kerrison Punch 5mm x 200mm @45°	36
S139	Long Laser Suction 6fg 20cm Flush Finish	21
S146	Nasal Aspirating Dissector	18
S147	Rounded Suction 12fg 9cm 4x1mm Fenestrations	22
S156	Rounded Long Suction 6fg 20cm	11
S157	Long Schuknecht 20g 12.5cm	18
S160	Fenestrated Suction 12fg 5cm	9

Reference	Description	Page
S161	Kerrison Punch Fine 1mm x 200mm for Delicate Use	36
S163	Sinus 12fg No Control Hole	17
S165	Laser Suction 6fg 9cm Flush Finish	21
S168	Fenestrated Rounded Suction 6fg 20cm @60°	12
S169	Fenestrated Rounded Suction 9fg 20cm @60°	12
S176	Suction 3fg 8cm	14
S177	Suction 5fg 8cm	14
S179	Thudichum Specula Size 5	38
S180	Tilley Nasal Forceps Short	38
S181	Tracheal Dilator	38
S182	Rounded Long Suction 9fg 27cm	11
S183	Rounded Long Suction 9fg 45cm	11
S184	Rounded Long Suction 12fg 55cm	11
S185	Rounded Long Suction 6fg 27cm	11
S187	Rounded Long Suction 9fg 35cm	11
S189	Teardrop Suction 9fg 9cm	13
S190	Teardrop Suction 6fg 9cm Tapered	13
S192	Teardrop Suction 9fg 13cm	13
S195	Diathermy Suction 12fg 19cm	24
S199	Suction 12fg 17cm	7
S205	Cushings Artery Forceps Curved, Half Serrated Jaws, 150mm	33
S206	Jefferson Retractor, Self-Retaining, Curved, Blunt Pointed Blades, 125mm	33
S207	Mollinsons Retractor, Self-Retaining, Curved, 4/4 Teeth, Adult, 165mm	33
S211	Hudson Brace, Standard Size with Universal Fitting, 265mm	34
S212	Gigli Saw 300mm and Handles, Pair, 55mm with Guide	34
S213	6fg Cannula Curved Set of 2	29
S214	6fg Cannula Straight	29
S221	Micro Ear Forceps 5mm	41
S228	Ball End Sinus 9fg	17
S229	Ball End Sinus 12fg	17
S230	Ball End Sinus Small Bend 9fg	17
S231	Ball End Sinus Small Bend 12fg	17
S236	Rounded Suction 12fg 5.5cm	8
S238	Thudichum Specula Size 4	38
S239	Thudichum Specula Size 3	38
S240	Sheffield Suction 9fg 15cm 90° Bend with Clip	20
S255	Iris Scissors Straight	35
S256	Iris Scissors Curved	35
S257	Fine Suture Forceps	35
S258	Long Suction 6fg 45cm	11
S261	Kerrison Punch Fine 2mm x 200mm for Delicate Use	36
S262	Kerrison Punch Fine 3mm x 200mm for Delicate Use	36
S263	Tischler-Morgan Biopsy Punch 200mm, Straight Tip	32
S284	Rounded Long Suction 12fg 35cm	11

Reference	Description	Page
S285	Fairley Laser Suction 12cm with Extended Angulated Tip	21
S301	Myringotomy Kit 1 Contains 11 items	39
S303	Jobson Horne Probe	40
S304	Ear Specula 4mm	42
S305	Ear Specula 5mm	42
S306	Ear Specula 6mm	42
S307	Cerumen Hook	40
S308	Rosen Insertter	41
S309	Ear Specula 3mm	42
S310	Myringotomy Knife	39
S318	Rounded Long Suction 12fg 20cm	11
S319	Rounded Long Suction 9fg 20cm	11
S320	Tischler-Morgan Biopsy Punch 200mm, Tip Down	32
S321	Soft Tip 6fg 9cm Suction	14
S340	Frazier Suction 15fg 9cm	7
S343	Fenestrated Rounded Frazier 12fg 13cm	10
S353	Pulse Lavage Kit includes Pulse Lavage and 2 tips	43
S367	Kerrison Punch 2mm x 200mm @90°	36
S380	Tilley Nasal Forceps Long	38
S419	Yasargil Retractor	34
S432	Insulated Suction 9fg 17cm 90°	23
S433	Insulated Suction 6fg 17cm 90°	23
S434	Micro Suction 18g 7cm	16
S435	Retractor Suction 12fg 12cm 60°	15
S436	Dissector Suction Flattened Tip 12fg 13cm	15
S438	Diathermy Suction 9fg 17cm 70°	25
S439	Teardrop Micro Suction 18g	13
S440	Poole Tip Suction 12fg 13cm Rounded 4x4x2mm Fenestrations	19
S444	Ball End Sinus Suction 9fg 60° Curve	17
S445	Rounded Suction 12fg 30cm 4x1mm Fenestrations 60°	22
S446	Curved Micro Suction 9fg 15cm 30° Bend 60° Curved Tip (Left)	16
S447	Curved Micro Suction 9fg 15cm 30° Bend 60° Curved Tip (Right)	16
S448	Diathermy Abbey Needle with Suction 6fg 10cm 40° Bend	26
S467	Ball End Suction Upward Curve 70° 9fg	17
S468	Ball End Suction Upward Curve 65° 12fg	17

CONTACT US

Eakin Surgical Ltd
Greypoint
Cardiff Business Park
Cardiff
CF14 5WF
United Kingdom

Telephone	+44 (0)2920 767 800
Facsimile	+44 (0)2920 767 801
Email	info@eakinsurgical.com
Website	www.eakinsurgical.com

Please send orders to orders@eakinsurgical.com

Please note products marked “*Made to Order” are not subject to standard delivery times and require a minimum order quantity. For more information, please contact info@eakinsurgical.com

These devices are intended for an individual patient, during a single procedure and should be discarded after use. Potential of cross contamination if re-used due to device not being sterile.

Eakin Surgical Ltd is registered in England & Wales. Company Registration No. 3939863. VAT No. GB806634432
Registered office: Eakin Surgical Ltd, Greypoint, Cardiff Business Park, Cardiff, CF14 5WF.

V16
09/2022